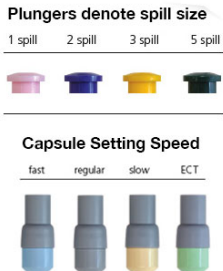


MATERIALS USED IN THE STUDENT CLINIC

	Dentine conditioner Brand: GC Composition: 10% polyacrylic acid (PAA) Use: 15 sec scrub, 10 sec wash, dry without desiccating	- Removes protein, oil and debris (smear layer) - Forms smear plugs to block dentine tubules - Prepares dentine for ion exchange adhesion
	Etchant Brand: ESPE Scotchbond Composition: 37% phosphoric acid Use for enamel: 30 sec scrub, 20 sec wash, 10 sec dry Use for dentine: 15 sec scrub, 20 sec rinse, 10 dry without desiccating	- Preferentially etches enamel - Forms tags of 10-30 µm - Increases the surface area for bonding - Partly demineralises dentine - Forms collagen "fibrils" - Prepares dentine for hybrid layer
	4th gen primer and adhesive Brand: ESPE Scotchbond Multi-Purpose Use: Apply primer to enamel and dentine, gently air blow, apply adhesive and allow 15 sec to penetrate hybrid layer, 20 sec light cure	- Provides a reliable bond to dentine - Prof Bryant's preferred system - Use after etching
	5th gen primer/adhesive Brand: ESPE Adper Single Bond Use: Shake bottle well, apply two consecutive coats to enamel and dentine, 5 sec gentle dry, 15 sec penetration time, 20 sec light cure	- Use after etching - Saves time vs. 4th gen system
	6th gen etchant/primer/adhesive Brand: ESPE Scotchbond Universal Use: Shake bottle well, 20 sec scrub, 5 sec gentle dry, 15 sec penetration time, 10 sec light cure	- One-step etchant, primer and adhesive - Saves time vs. 5th gen system - Preferred method at CSU is to use in conjunction with selective enamel etching (modified 6th gen)

	RMGI liner Brand: ESPE Vitrebond Use: Thoroughly mix one spoonful of powder with one drop of liquid and use immediately	• Seals exposed dentine tubules • Radiopaque • Thick consistency, can replace lost dentine • Favoured beneath CR • Bulk reduces polymerisation contraction • Protects pulp from LMW resin components • 2-3 minute working time
	RMGI liner Brand: GC Fuji Bond Use: Thoroughly mix one spoonful of powder with two drops of liquid and use immediately	• Seals exposed dentine tubules • Radiopaque • Thin consistency, allows bulk of restorative material • Favoured beneath amalgam • Protects pulp from LMW resin components • 2-3 minute working time
Fuji IX (GIC) • Fluoride release • Ion-exchange bonding to dentine Indications • Core material • Non-aesthetic posterior Class I, II and V restorations • Base layer in Class II sandwich restorations Components • 95% Glass base • 5% PAA		
	GIC Brand: GC Fuji IX (grey button) Use: • 10 sec mixing time • 2 minutes working time • 4 ½ minutes setting time	• Standard formula
	GIC Brand: GC Fuji IX Fast (white button) Use: • 10 sec mixing time • 1 minute working time • 3 minutes setting time	• Fast set • Challenging patients
	GIC Brand: GC Fuji IX Extra (yellow button) Use: • 10 sec mixing time • 1 minute working time • 2 ½ minutes setting time	• Extra fast set • Challenging patient • 6x higher fluoride than Fuji IX Fast • High compressive strength • Extra translucency • Geriatrics and paed
	RMGI Brand: GC Fuji II LC (green button) Use: • 10 sec mixing time • 3 minutes working time • 20 sec light cure	• Dual cure • Class III • Class V • Deciduous teeth • Liner • Base • Core build-up

	RMGI Brand: GC Fuji VII (pink button) Use: • 10 sec mixing time • 2 ½ minutes working time • 4 minutes setting time (GI component)	• Dual cure • Enhanced with CPP-ACP • Fissure sealants in deciduous teeth (pink) • Fissure sealants in permanent teeth (white) • Control of hypersensitivity • Temporary filling material • 6x higher fluoride than other sealants
	Composite resin Brand: GC Gradia Direct Use: Highly aesthetic restorations	• Microfilled • 'Sticky' • Difficult handling characteristics • Good translucency
	Composite resin Brand: ESPE Filtek Supreme XTE Use: Restorations with complex increments, including aesthetic regions	• Nanofilled • Easier handling characteristics • Average translucency
 Plungers denote spill size 1 spill 2 spill 3 spill 5 spill Capsule Setting Speed fast regular slow ECT	Amalgam Brand: SDI Permite C Use: Large, non-aesthetic restorations	• Spherical and lathe cut particles • High Cu admixed alloy (Cu, Ag, Sn) • High compressive strength • Corrosion products seal margins • Only 2-spill regular and ECT available at CSU
	Topical anaesthetic Brand: Astrazeneca Xylocaine Use: Topical anaesthesia prior to injection	• 100 mg/g lignocaine • Onset: 3-5 minutes • Duration: 20 minutes • Maximum dose: - o Adult: 5 mg/kg - o Child: 3 mg/kg
	Topical anaesthetic Brand: Dentsply Oraqix Use: Topical anaesthesia of periodontal pockets	• 25 mg/g lignocaine • 25 mg/g prilocaine • Changes from gel to liquid at body temperature • Onset: 30 seconds • Duration: 20 minutes • Maximum dose: 5 cartridges/session
	Injectable anaesthetic solution Brand: Septodont Lignospan Special Use: Local anaesthesia, with vasoconstriction	• 20 mg/mL lignocaine • 12.5 µg/mL adrenaline • Onset: 1-4 minutes • Duration: 1-2 hours • Maximum dose: 7 mg/kg

	Injectable anaesthetic solution Brand: Septodont Scandonest Plain Use: Local anaesthesia, without vasoconstriction - Metabisulphite allergy - Hypertension	3 mg/mL mepivacaine - Onset: 4-5 minutes - Duration: 2-3 hours - Maximum dose: - Adults: 6.6 mg/kg 6-14 years: 2.7 mL 3-6 years: 1.8 mL
	CaOH pulp capping material Brand: Kerr Life Use: Mix catalyst and base in equal amounts. Use immediately.	- Indirect or direct pulp capping - Anti-microbial due to high pH (11-13) - Hard set - Must be sealed with RMGI - May stimulate reparative dentine
	Steroid-based pulpal agent Brand: Sigma Ledermix Paste Use: Temporary obturation of root canals	10 mg/g triamcinolone acetonide - Corticosteroid 30 mg/g demeclocycline HCl - Broad spectrum antibiotic - Water soluble
	Steroid-based pulp capping material Brand: Sigma Ledermix Cement Use: Mix powder and liquid (ratio not critical) and use immediately	- Powder 6.7 mg/g triamcinolone acetonide 20 mg/g demeclocycline HCl 47.2% zinc oxide 33.4% CaOH - Liquid 85% eugenol - Hard set
	Steroid-based pulpal agent Brand: ADM Odontopaste Use: Temporary obturation of root canals	10 mg/g triamcinolone acetonide - Corticosteroid 50 mg/g clindamycin HCl - Broad spectrum antibiotic
	Steroid-based pulp capping material Brand: ADM Odontocem Use: Mix liquid and powder in approx. 3:1 ratio and use immediately	2 mg/g triamcinolone acetonide - Biocompatible calcium silicate cement - Resorbable
	Haemostatic agent Brand: Ultradent Viscostat Use: Agitate in areas of bleeding, wait 30 seconds and rinse coagulation. Repeat if necessary to achieve haemostasis	200 mg/g ferric sulphate - Potent haemostatic agent - Coagulation can stain
	Fluoride varnish Brand: Colgate Duraphat Use: Apply directly to demineralised areas	- Fluoride content = 22,600 ppm - Don't eat or drink for ½ hour - Varnish will drop off at some point - Safe to ingest - Only for 10+ years of age